

# FDA Submission

**Your Name:** Deep Learning Researcher

**Name of your Device:** PneumaTron-X

## Algorithm Description

### 1. General Information

**Intended Use Statement:**

This device is a deep neural network algorithm which is intended for classification of presence or absence of pneumonia in the age group of 1-100 years irrespective of any gender or racial bias.

**Indications for Use:**

To assist radiologist in classification of presence or absence of pneumonia

**Device Limitations:**


System requirement at the least is a CPU with 8GB RAM and 50 GB memory. Though having a GPU can improve the operational efficiency.

**Clinical Impact of Performance:**

This device was designed to achieve a F1 score higher than 0.3, leading to a higher Recall rate. This will lead to a higher False positive rate, which means it can classify a prediction as Positive despite it being negative for Pneumonia. But this will assist radiologists get better at their X-ray readings with a good accuracy for positive tests.

## 2. Algorithm Design and Function

### Algorithm Flowchart

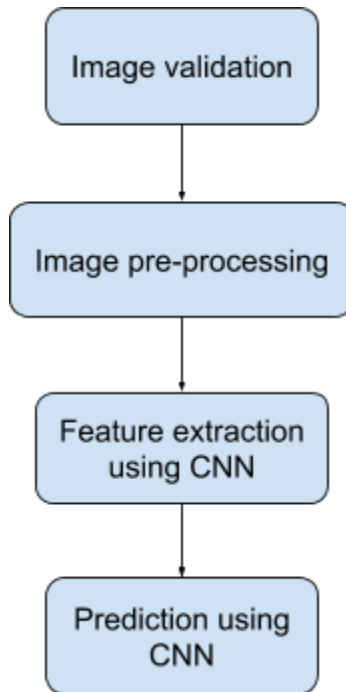

### DICOM Checking Steps:

DICOM checking involves three tests

- Modality test to ensure image is of Digital X-ray type
- Patient position to ensure image has Anteroposterior(AP) or Posterior-Anterior(PA) position of patient
- Patient body to ensure image is a Chest X-ray

Once image passes all three DICOM checking steps, it is feed to the device

### Preprocessing Steps:

First part of the device is pre-processing which ensures images fed to the deep neural network algorithm are processed well to ensure best algorithmic performance. This involves image resizing to 224 by 224 pixels, and normalized from 0,255 to 0,1

## **CNN Architecture:**

The CNN architecture is based on VGG16 base model, and has following layers on the top of [VGG16 model](#):

- Flatten
- Dropout with 50% probability
- Dense, 1024 units with ReLU activation
- Dropout Dropout with 50% probability
- Dense, 512 units with ReLU activation
- Dropout with 50% probability
- Dense, 256 units with ReLU activation
- Dense, 1 unit with Sigmoid activation

### 3. Algorithm Training

Parameters:

**Types of augmentation used during training:**

Data augmentation was only applied to the training set images, which included:

- A horizontal flip
- Random height shift of  $\pm 10\%$  of image height
- Random width shift of  $\pm 10\%$  of image width
- Random rotation transformation of  $\pm 20$  degrees
- Random shear transformation of  $\pm 10$  degrees,
- Random zooming of  $\pm 10\%$

**Batch size:** 32 data points per batch

**Optimizer learning rate:** Adam optimizer with learning rate of 0.001 was used while training

**Layers of pre-existing architecture that were frozen:**

All the layers except last output layer of the base VGG16 model was frozen, 17 in total


**Layers of pre-existing architecture that were fine-tuned:**

Last output layer of base VGG16 model was fine tuned to flatten it for new layer addition


**Layers added to pre-existing architecture:**

Already described in CNN architecture

## Algorithm training performance visualization

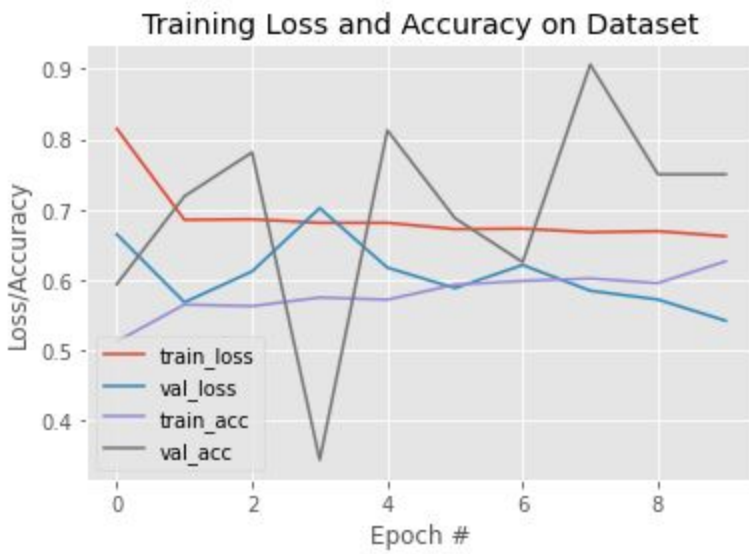

## P-R curve

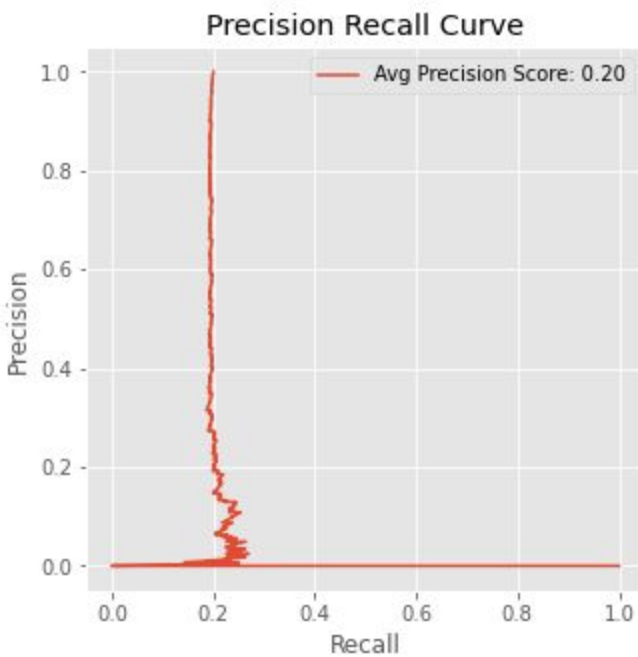


## Final Threshold and Explanation:

To get a F1 score higher than 0.3, a threshold of 0.349681 was selected.

## 4. Databases

National Institute of Health(NIH) Chest X-rays dataset was used for this project. This dataset is hosted on the [Kaggle website](#). There are 112,120 X-ray images with disease labels from 30,805 unique patients in this dataset. The labels include 14 common thoracic pathologies:

- Atelectasis
- Consolidation
- Infiltration
- Pneumothorax
- Edema
- Emphysema
- Fibrosis
- Effusion
- Pneumonia
- Pleural thickening
- Cardiomegaly
- Nodule
- Mass
- Hernia

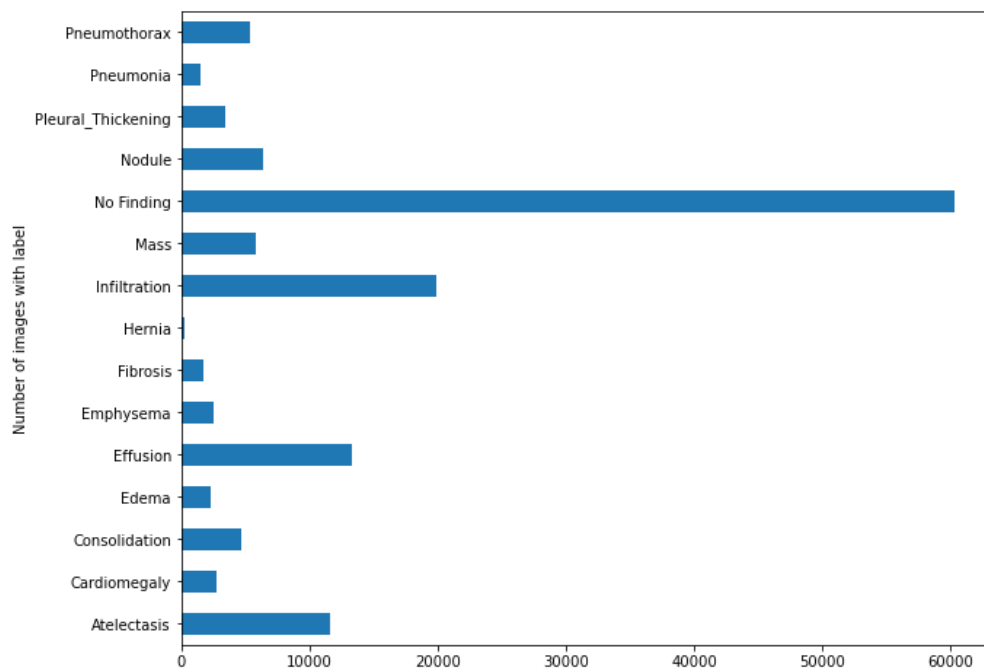



Disease wise distribution of the dataset can be seen in the above figure.

## Description of Training Dataset

Training dataset had 1145 number of pneumonia cases and 88551 number of non-pneumonia cases. Percentage of pneumonia cases in the train set was 1.28% which was balanced to 50% by dropping non-pneumonia cases. Total representation of this updated train set:

Male: 1325

Female: 965

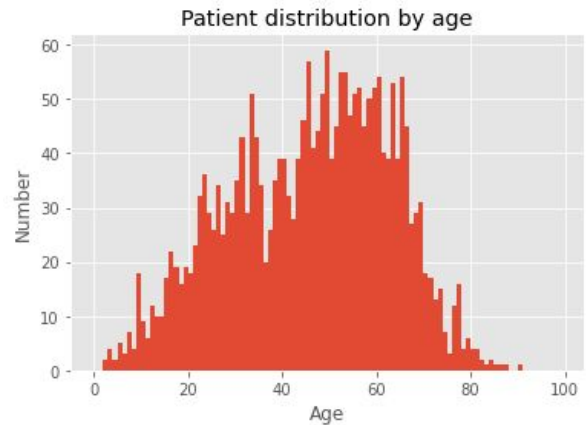
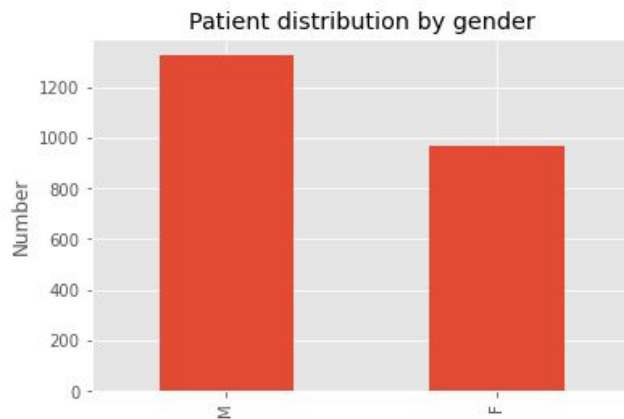


## Description of Validation Dataset

Validation dataset had 286 number of pneumonia cases and 22138 number of non-pneumonia cases. Percentage of pneumonia cases in the validation set was 1.28% which was balanced to 20% to simulate the real world scenario by dropping non-pneumonia cases. Total representation of this updated train set:

Male: 810

Female: 620

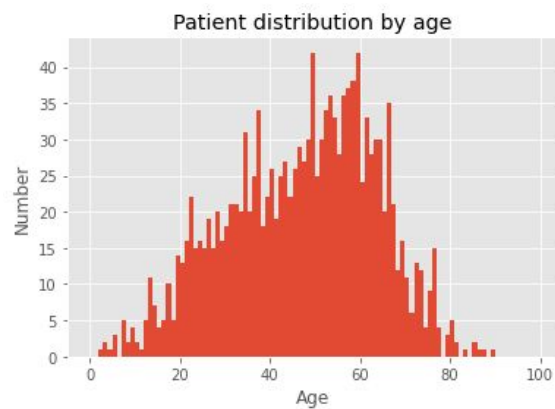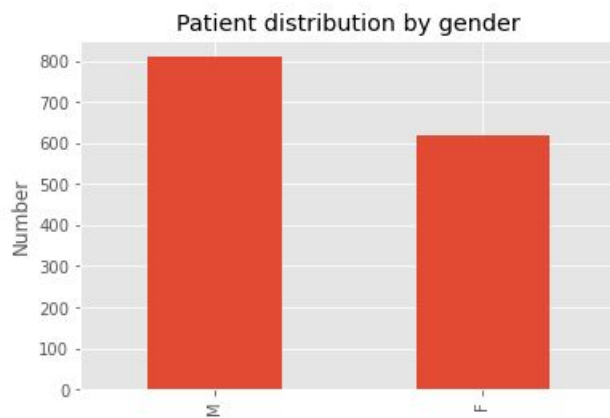

## 5. Ground Truth

Ground truth for the disease labels was created using Natural Language Processing (NLP) to mine the associated radiological reports. Since the labels are NLP-extracted, there could be some erroneous labels but the NLP labeling accuracy is estimated to be more than 90%

## 6. FDA Validation Plan

### **Patient Population Description for FDA Validation Dataset:**

For FDA validation set, the partner hospital or institute should create a dataset with both male and female cases in almost equal number, spread across age groups of 1 to 100 years. Dataset can contain patients with or without prior diagnosis of any of 14 common thoracic diseases.

### **Ground Truth Acquisition Methodology:**

The ground truth should be acquired in AP or PA position of the chest part of a patient and should be of a Digital X-ray type.

### **Algorithm Performance Standard:**

Algorithm performance can be measured in terms of F1 score. Recall and Precision can also be used to measure standard performance, but F1 score will give an accurate representation of the performance.